




FaceFit™ NIV masks

Lightweight and comfortable NIV (Non-invasive Ventilation) masks for conventional ventilation or CPAP and bilevel applications

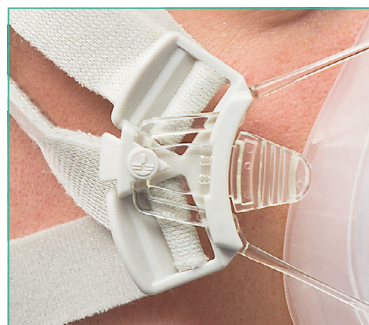

Critical Care ▪ Patient Interfaces



FaceFit™ NIV masks

FaceFit is a comfortable single use adult hospital mask for use in non-invasive ventilation.

Two options are available; one option is for use with conventional dual limb breathing systems, and the second for bilevel and CPAP applications where an anti-asphyxiation valve and CO₂ leak port is required.

FaceFit masks are available in small, medium and large sizes to fit a wide range of patients. A size guide is included to help with making the right choice prior to opening the packaging.

Adjustable rear head pad

Allows for optimal fit for a wide range of patients

Lightweight

Improving patient tolerance

Easily adjustable harness and quick release mechanism

Allowing for fast removal of the mask from the patient and optimum fit

Video available at www.intersurgical.com



Anti-asphyxiation valve (FaceFit ported only)

Prevents rebreathing in the case of inadequate ventilation

CO₂ Leak port (FaceFit ported only)

For effective CPAP or bilevel ventilation

Low dead space

For efficient respiratory support

6mm stem connector with retainable cap

For supplementary oxygen or monitoring

Soft silicone seal

Ensuring optimum performance

illustrated

FaceFit mask without anti-asphyxiation valve and without CO₂ leak port

For use with dual limb breathing systems

Code	Description	Box Qty.
2255000	FaceFit, NIV mask, small adult	12
2256000	FaceFit, NIV mask, medium adult	12
2257000	FaceFit, NIV mask, large adult	12

FaceFit ported with anti-asphyxiation valve and CO₂ leak port

For use with single limb breathing systems that do not have a CO₂ leak port

Code	Description	Box Qty.
2250000	FaceFit ported, NIV mask with anti-asphyxiation valve and CO ₂ port, small adult	12
2251000	FaceFit ported, NIV mask with anti-asphyxiation valve and CO ₂ port, medium adult	12
2252000	FaceFit ported, NIV mask with anti-asphyxiation valve and CO ₂ port, large adult	12

Make an enquiry

IS1.15 • Issue 4 03.17



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
T: +44 (0)118 9656 300 F: +44 (0)118 9656 356 info@intersurgical.com www.intersurgical.com

UK • Ireland • France • Germany • Spain • Portugal
Italy • Benelux • Sweden • Denmark • Lithuania • Russia
Czech Republic • Turkey • South Africa • China • Japan
Taiwan • Philippines • USA • Canada • Colombia • Australia



Please think before you print
Save energy and paper.
If you must print this e-mail
please print it double sided

Intersurgical Ltd is certified to ISO 14001, ISO 9001 and ISO 13485